


Setup Scenario and Measurement

(August 09, 2016, MHK 628, 3:00 pm. to 7:40 pm.)

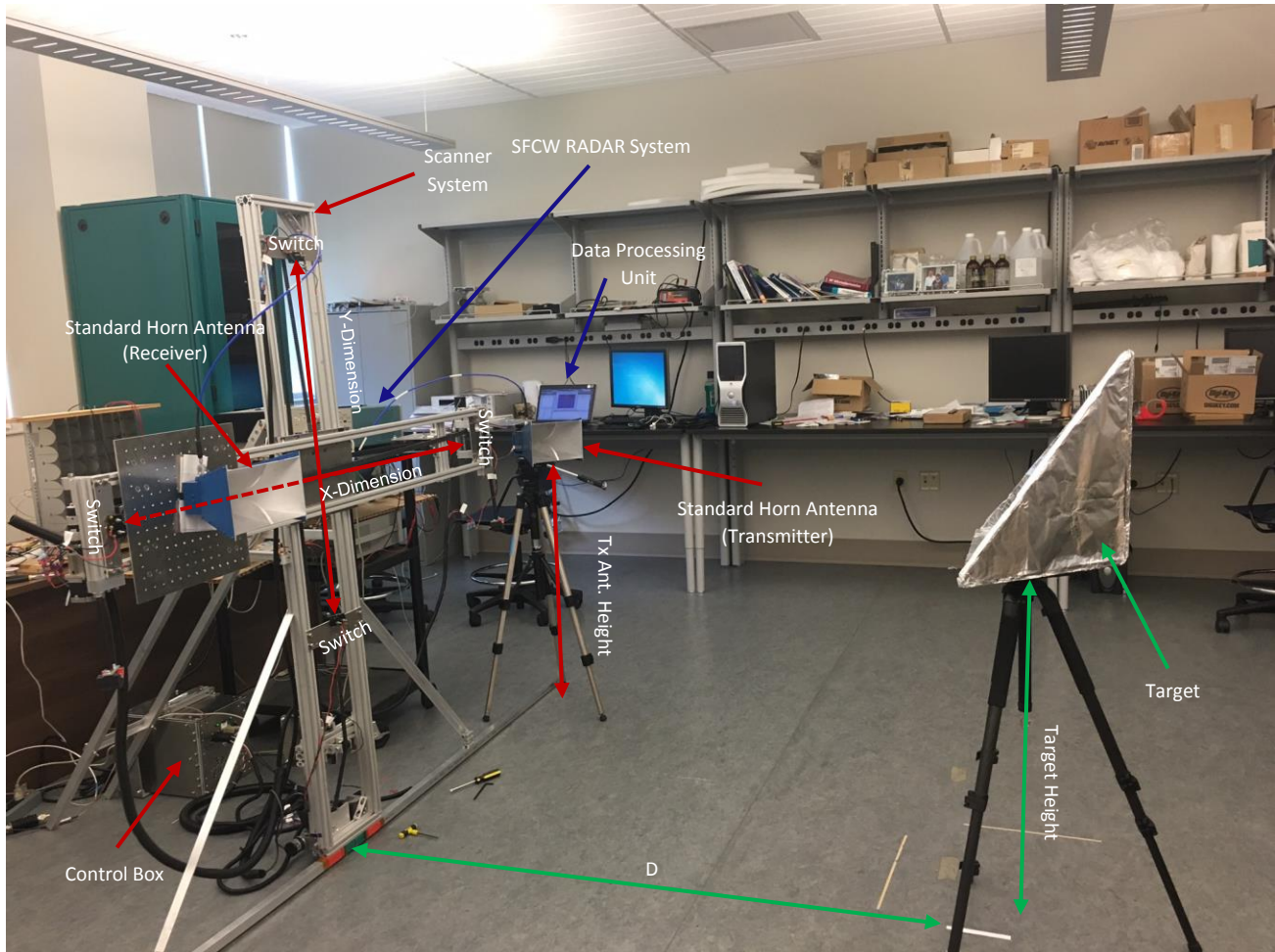

Control Box-Software Setup:

Scan Area:

X-Dimensions (cm): 122

X-Step: 5 cm

Dwell Time: 10 sec.

Y-Dimensions (cm): 113

Y-Step: 5 cm

X-Dimension & Y-Dimension are the distances between 2 switches

D: 160 cm

Target Height: 92 cm

Tx Ant. Height: 105.5 cm

Rx Ant. Height @ Origin Point: 84.5 cm


Distance between Tx Ant. and Target: 180 cm

Distance between Rx Ant. @ Origin Point and Target: 160 cm

Minimum Distance between Tx and Rx Antennas: 32 cm

Distance between Tx and Rx Antenna @ Origin Point: 168 cm

Scan Pattern (Stop Points) (Front View) of Receiving Antenna (Standard Horn Ant.)

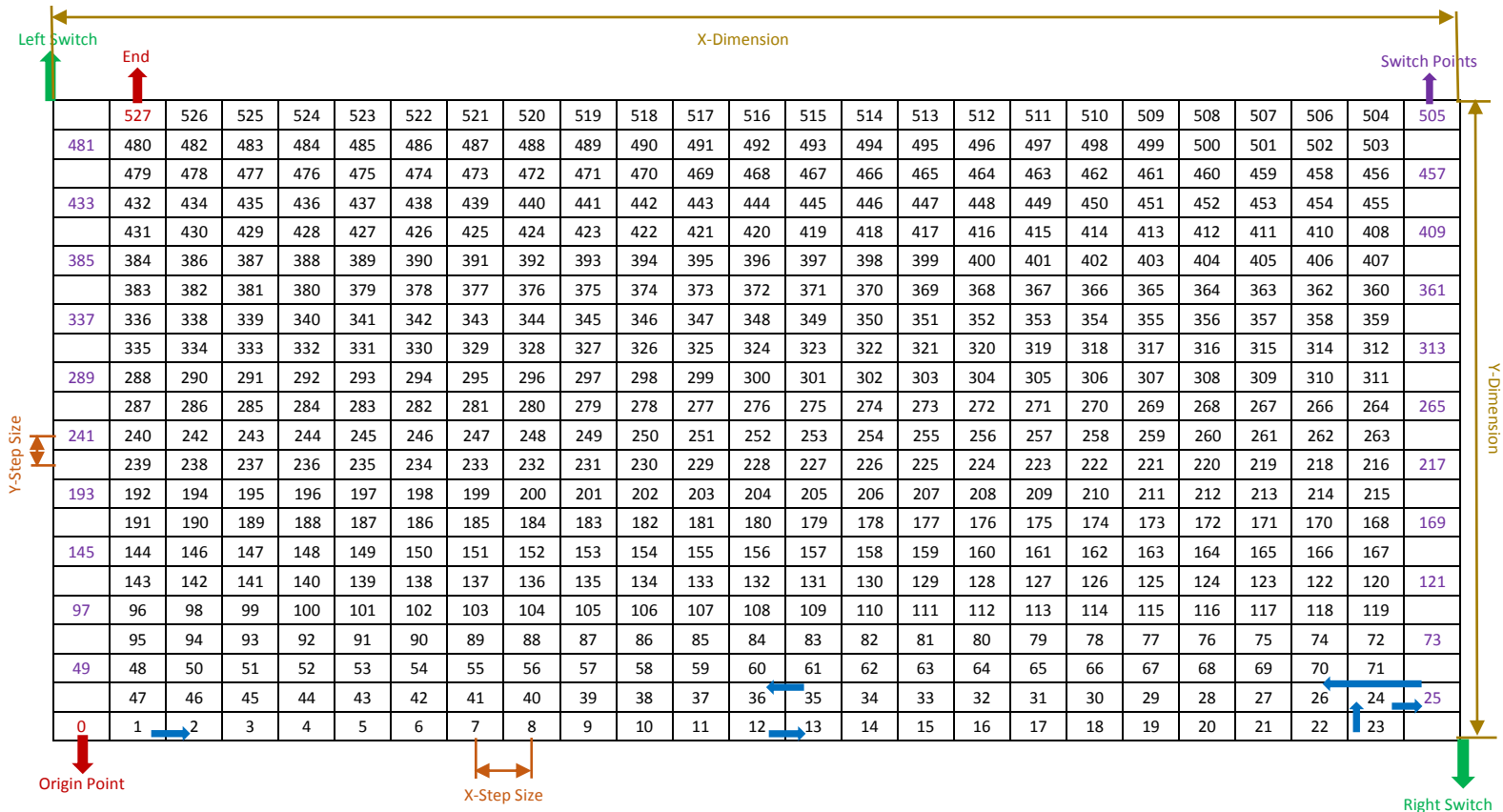

Dwell Time (Stop Time) @ Each Stop Point: 14sec.

Dwell Time (Stop Time) @ Each Switch Point: 28 Sec.

LabView Data-Collecting Time: 1 sec. (We collect data at each stop point for 1 seconds using Labview.)

Real Scan Area:

X-Dimensions (cm): 158.5

X-Step: 5 cm

Dwell Time: 14 sec.

Y-Dimensions (cm): 96

Y-Step: 3 cm

Plate Size: 35.5 cm × 35.5 cm

X-Dimension & Y-Dimension are the distances between 2 switches

Notes:

- Every time we close and open the software (to send commands to control box and scanner system), we need to push Startup System Button first and wait until the plate comes to origin point (Then we can initialize setup and push Start Scan button).
- While the system is scanning, if we turn the system off using the ON/OFF switch on control box (but not close the controlling software) and turn it on again, it will start scanning from the last position (but the software must be open).
- When the system finishes the scanning, we can push Go To Origin button and the plate will return to the Origin Point from the End Position (The system should be on and the software should be open).